




stella

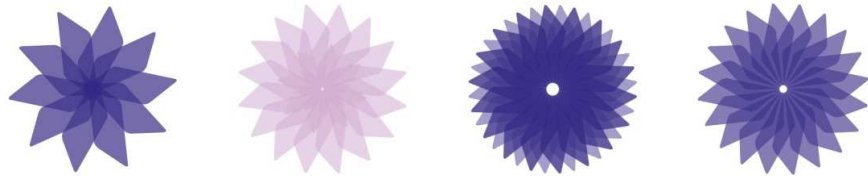

Stella Suite Bluetooth Installation User Guide

Version 1.4
© Copyright Reserved 2011 Tristel Solutions Ltd



Stella Suite

Bluetooth Installation User Guide

Table of Contents

1.	Setting up Bluetooth	3
2.	Download via Bluetooth	6
3.	Product Support Service	7



1. Setting up Bluetooth

Connecting Stella IQ to Stella Suite via Bluetooth


By initiating a Bluetooth connection on Stella IQ, events can be downloaded to Stella Suite software without the use of a USB cable.

1.1 From your computer select the start menu/ control panel/ devices and printers, then select “Add a device”.

1.2 Set Stella IQ to download mode, by pressing the “ON” button twice. Use the “ON” button to toggle down to the “Download” option. Stella IQ will now automatically start a Bluetooth establishment process.




1.3 “Add a Device” window on your computer will show a Stella ICON followed by a two or three numbers. This is the serial number of the Stella IQ unit that you are establishing a Bluetooth connection with.

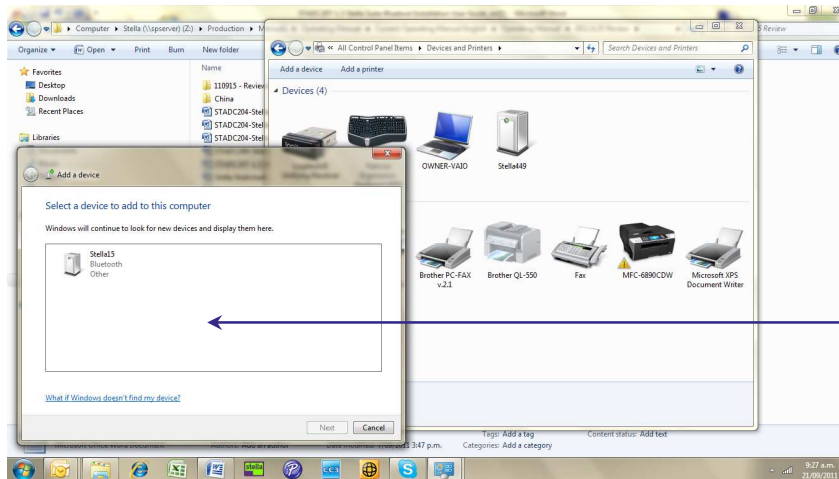


Stella IQ device added



- 1.4 Double click on the Stella device icon and the “select pairing option” window will open. Choose the second option “enter the devices pairing number”, enter 1234 then select next.

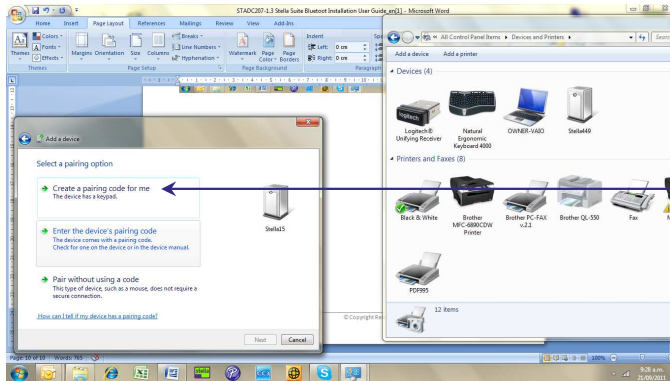


Select 'enter the device pairing numbers' option

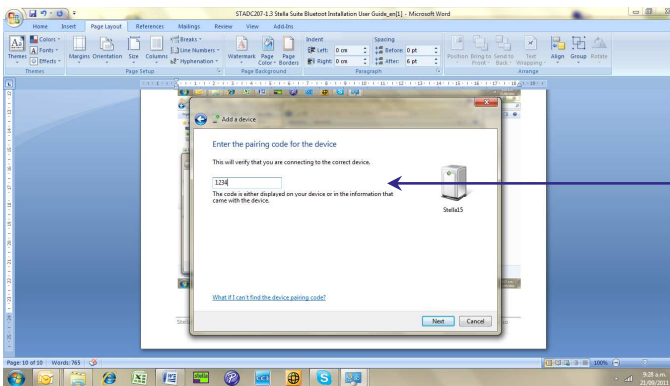

Entering pairing code 1234


- 1.5 From the control panel device and managers window, right click on the Stella device ICON and select properties. This will then provide a COM Port number such as COM 3. Take note of this number.

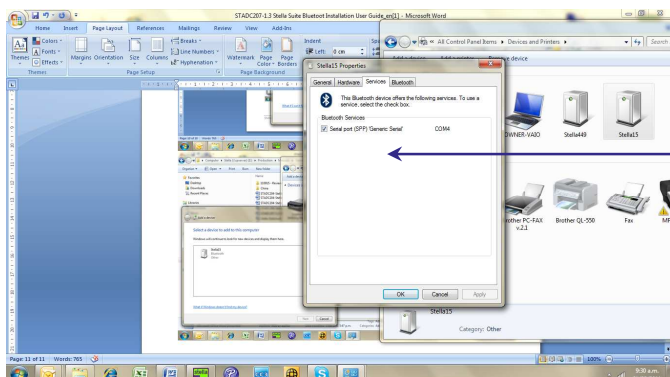


Take note of COM Port number



- 1.6 Close the Control Panel windows and open Stella Suite. Go to the “File” drop down menu and select ‘configure’ then the Bluetooth tab. Tick the ‘Use Bluetooth’ option and select the COM port number noted when installing the device. Select OK.

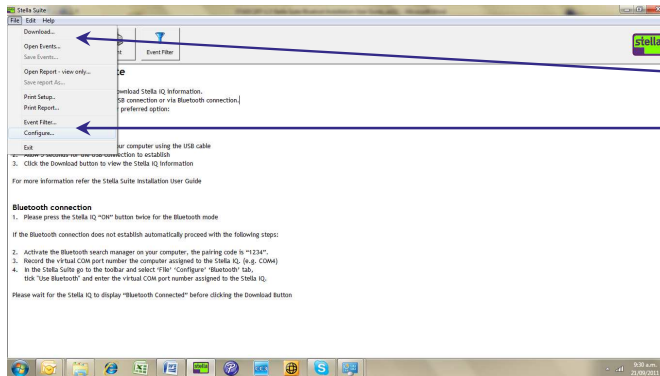

“File” menu
“Configure” window

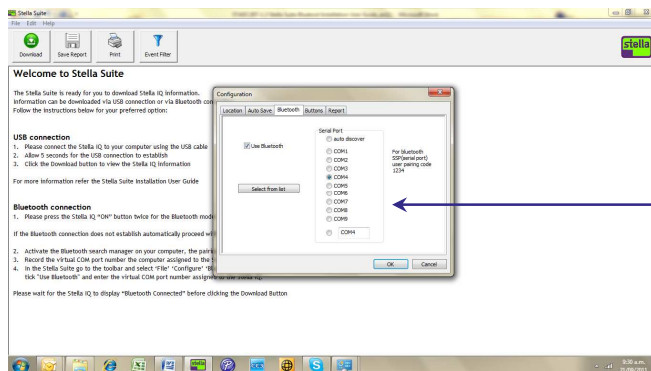

Select correct COM Port number
recorded earlier

- 1.7 The following graphic is displayed while the Bluetooth connection is being established.

Bluetooth is now active and ready to download.

- 1.8 To download press “the Download” button. Stella will now download an event file and data log report.

⚠ Bluetooth setup must be established for each Stella IQ unit.






2. Download via Bluetooth


Once Bluetooth set up is achieved, downloading events can be utilised by the following steps.

2.1 Open Stella Suite on your computer

2.2 Turn the Stella IQ “ON” button twice and toggle using the “ON” button to select the download option. When the following Icon appears, you are ready to download.

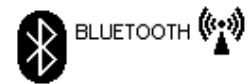

2.3 Select Download on Stella Suite. The Suite will now start a download process window

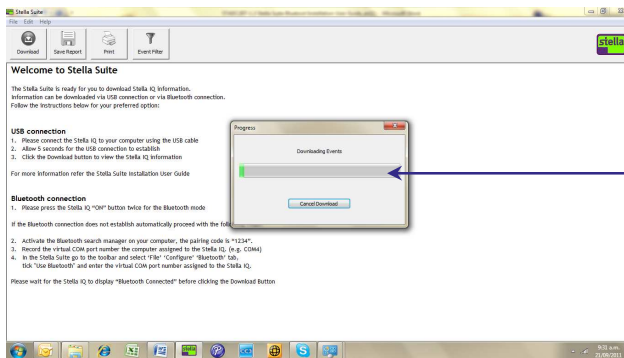

Download in progress

2.4 Stella IQ will now display the following download graphic.



For further information on Stella Suite and its functions, please refer to the Stella Suite Installation and User guide.



3. Product Support Service

If you have questions regarding the operation of this product please read this installation and user guide carefully before making further enquiries.

For more information on this product or other products please contact your local Stella Sales Agent. Contact details for Stella Sales Agents can be found in the Contacts section of the Operating Manual.

Tristel Limited Offices

New Zealand

Tristel New Zealand Limited
23 Birch Avenue,
Judea, Tauranga
New Zealand
T: +64 7 5771560
F: +64 7 5771567
E: info@stella-performance.com
W: www.stella-performance.com

China

Shanghai Stella Medical Equipment Co. Ltd
Rm.602, No.26, Lane 168, Da Du River Road
Putou District
Shanghai 200062
China
P: 86 21 3250 9758
M: 86 136 3250 9757
E: info@stella-performance.com
W: www.stella-performance.com

United Kingdom and Rest of the World

Tristel Solutions Limited
Unit 1B, Lynx Business Park
Fordham Road
Snailwell, Cambridgeshire
CB8 7NY
UK
T: +44 (0)1638 722716
F: +44 (0)1638 722725
E: info@stella-performance.com
W: www.stella-performance.com



Manual Version

Date	Version	Author	Changes	Check	Approve
21/09/11	1.4	JBEALE	Complete revision following Windows 7 installation process	JB	WS
07/04/10	1.3	BKENT	Removed “5” from Stella5 IQ	B Kent	JB
05/05/09	1.2	BKENT	Revise Contact details, Front cover page, Meta Normal-Roman Font, Removed reference to Vista, amended graphics	JB	AH
04/02/09	1.1	DRODGES	Original Version	B Kent	AW

STADC207